

PELVI.LOC® PL-3DS for moulded Seat

PELVI.LOC® PL-3D Pelvis Belts

HMV.Nr.: 26.99.02.2000




Version

The pelvic pads are 3 - dimensional shape, for optimum fixation of the pelvis without pressure on the spine anterior superior. Formed shape with ultra thin end at the moulded seat fixing point to avoid lateral pressure marks at the pelvis.

Prevents lateral pressure points on the pelvis. Nylon strap for passing through the moulded seat. The ratchet will run in addition to a strap to prevent injuries.

Performance


Seat shell with **PL-3DS-S** pelvic retraction belt.

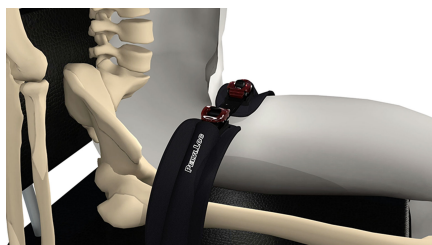

Amendment

For mounting on the moulded seat or on the wheelchair.

Article number: **PB-BL**

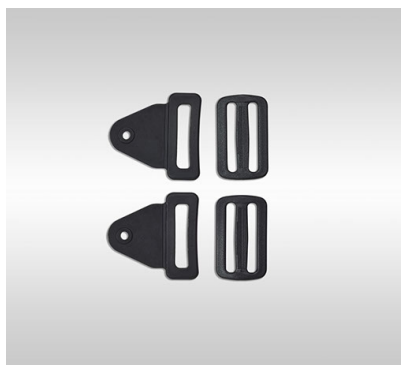

Scope of application

For custom made seats/ moulded seats.

Special features



Facet-like thinning of the padding material in the transition from the pad to the fastening strap. Prevents lateral pressure points on the pelvis.

Article numbers

PL-3DS-XXXX for SB 8-16 cm Bambino
PL-3DS-XXS for SB 12-20 cm
PL-3DS-XS for SB 16-24 cm
PL-3DS-S for SB 20-32 cm
PL-3DS-M for SB 30-40 cm
PL-3DS-L for SB 36-46 cm
PL-3DS-XL for SB 44-52 cm

Technical Data

	XXXX Bambino	XXS	Xs	S	M	L	XL
Seat width (cm)	8 - 16	12 - 20	16 - 24	20 - 32	30 - 40	36 - 46	44 - 52
Padding size (cm)	11 x 6	14 x 6	18 x 7	21 x 8	25,5 x 8	26,5 x 8,5	28,5 x 9
Weight (g)	220	220	220	220	280	280	280
Predetermined breaking point (at approx. kg)	40	40	60	60	80	80	80

Ratchet size XXXS - XXS	1,4 x 0,2 x 14,0 cm	<p>Product is latex-free</p> 	
Ratchet size XS - S	1,9 x 0,2 x 17,0 cm		
Ratchet size M - L - XL	2,4 x 0,2 x 17,0 cm		
Buckle size XXXS	B-14-STAR-BL		
Buckle size XXS	B-14-GO-BL		
Buckle size XS - S	B-19-ERGO-BL		
Buckle size M - L - XL	B-24-ERGO-BL		
Hole-strap			